




designed for scientists

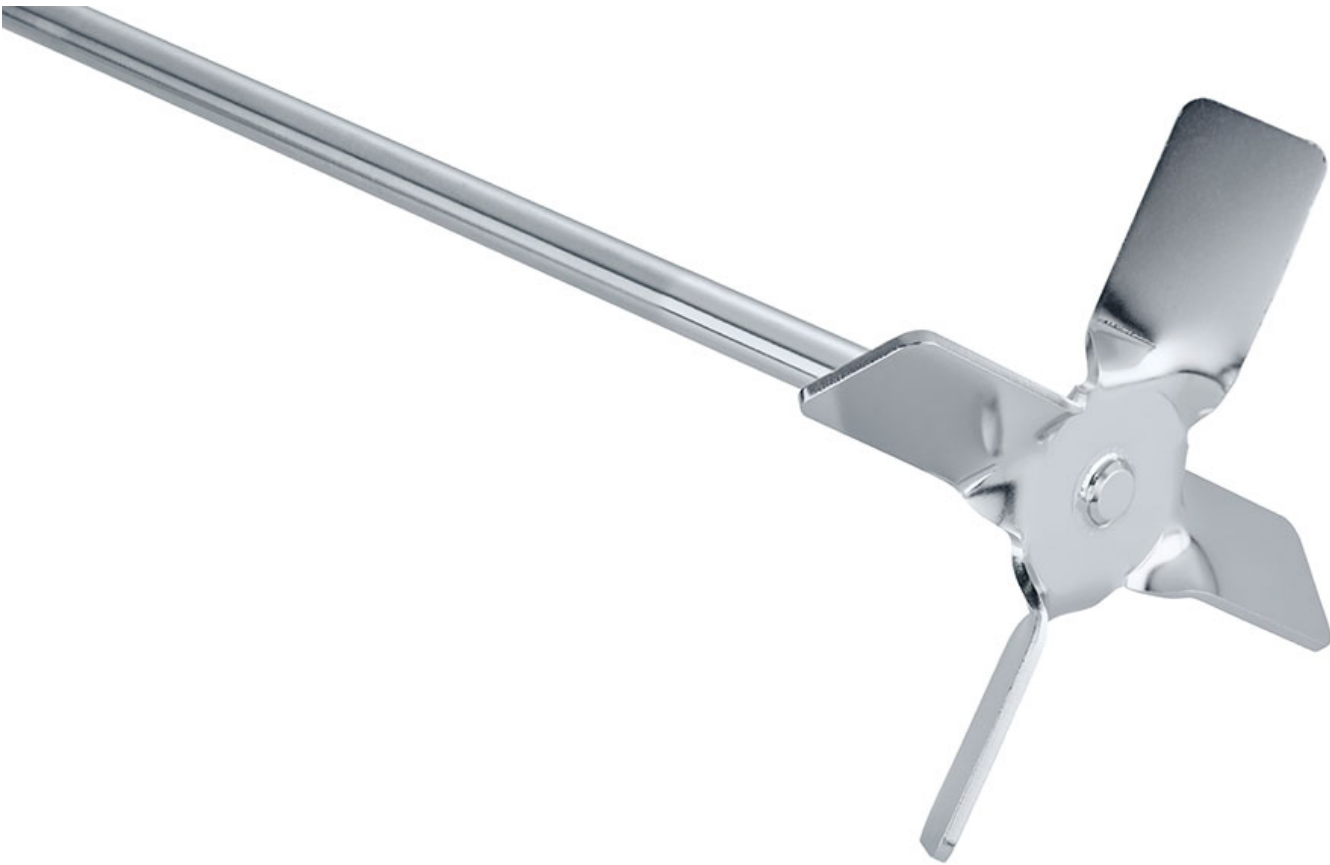

## R 2302 Propeller stirrer, 4-bladed

/// Data Sheet

Standard stirring element. For drawing the material to be mixed from the top to the bottom. Local shearing forces. Generates axial flow in the vessel. Used at medium to high speeds.

Stirrer Ø: 150 mm

Shaft Ø: 13 mm

Shaft length: 800 mm



designed for scientists

### Technical Data

Stirrer diameter [mm]	150
Shaft diameter [mm]	13
Shaft length [mm]	800
Speed max. [rpm]	600
Material in contact with medium	AISI 316L / 1.4571
Weight [kg]	1.024

