




designed for scientists

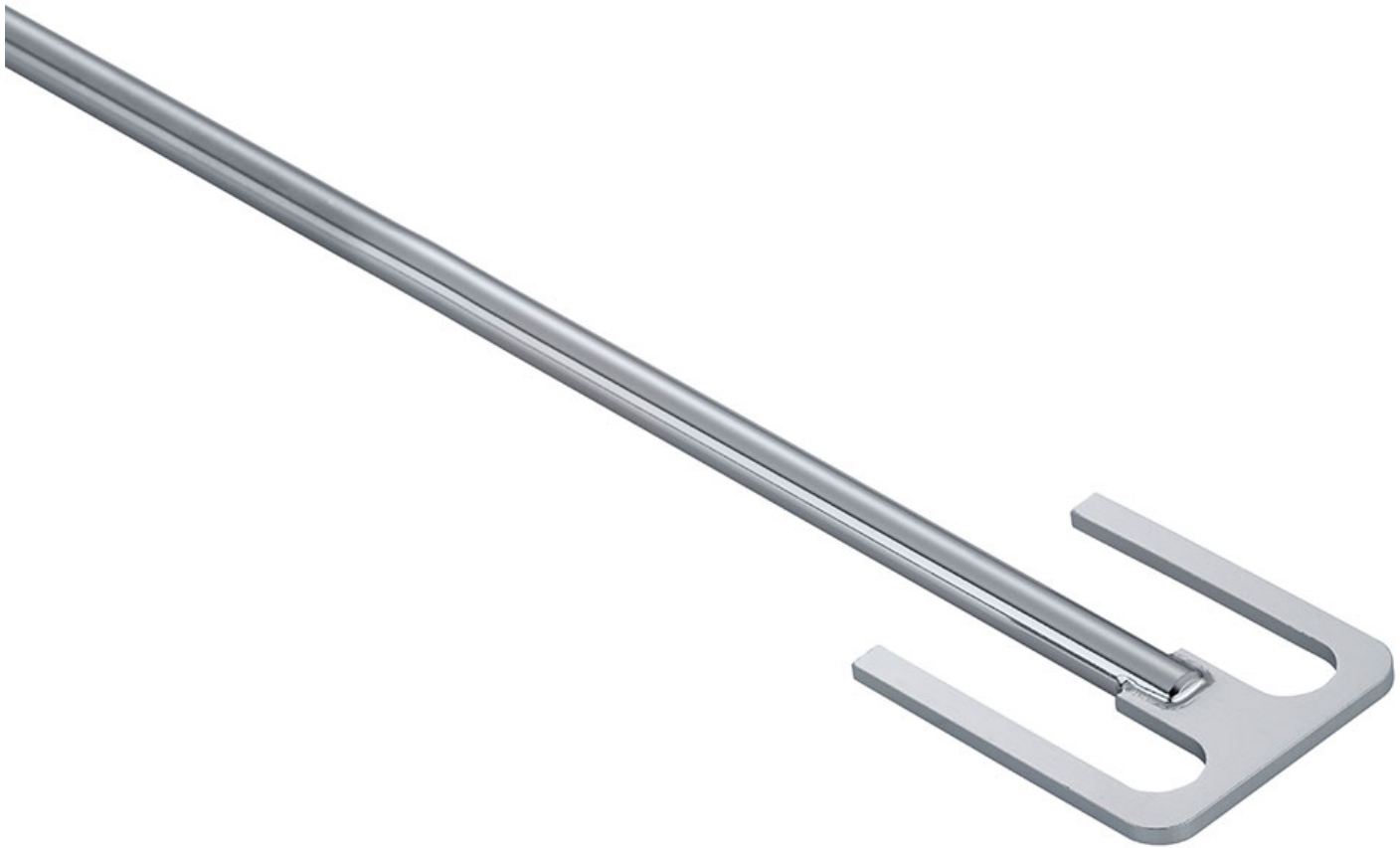

R 1330 Anchor stirrer

/// Data Sheet

Tangential flow, high shearing rate at edges, minimum deposits on the vessel wall. Used at low speeds. Polymer reactions, even distribution of high mineral contents in liquids. The ideal stirrer for medium to highly viscous fluids.

Stirrer Ø: 45 mm

Shaft Ø: 8 mm

Shaft length: 350 mm



designed for scientists

Technical Data

Stirrer diameter [mm]	45
Shaft diameter [mm]	8
Shaft length [mm]	350
Speed max. [rpm]	1000

