




designed for scientists

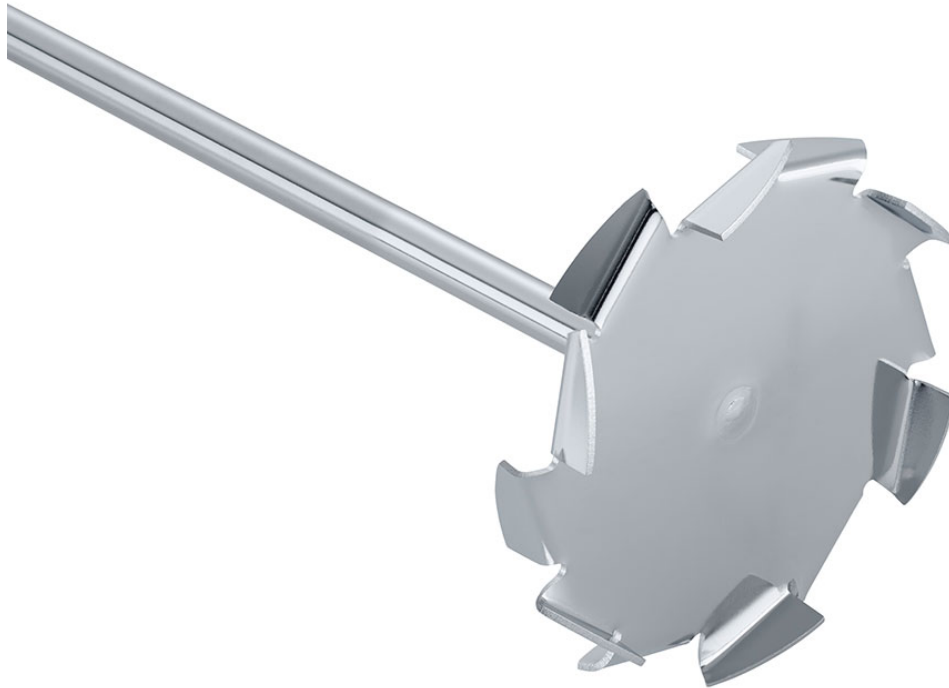

R 1302 Dissolver stirrer

/// Data Sheet

Radial flow, for drawing the material to be mixed from the top and the bottom. High turbulence, high shearing forces. Particle reduction. Used at medium to high speeds.

Stirrer Ø: 100 mm

Shaft Ø: 10 mm

Shaft length: 350 mm



designed for scientists

Technical Data

Stirrer diameter [mm]	100
Shaft diameter [mm]	10
Shaft length [mm]	350
Speed max. [rpm]	1000
Material in contact with medium	AISI 316L / 1.4571
Weight [kg]	0.312

