

# IKA

designed for scientists

EN




SMARTER WEIGHING | IKA [scale]

## IKA [scale]

/// The first stand scale of its kind

A stand with an integrated scale and data interface: only available from IKA. The IKA [scale] provides increased comfort and safety for your laboratory applications. For example, your stirrer can hang on the stable support rod. The quantities of samples to be added are weighed immediately and continuously documented during the stirring process.

**Patented:** The IKA [scale] is unique in the market!

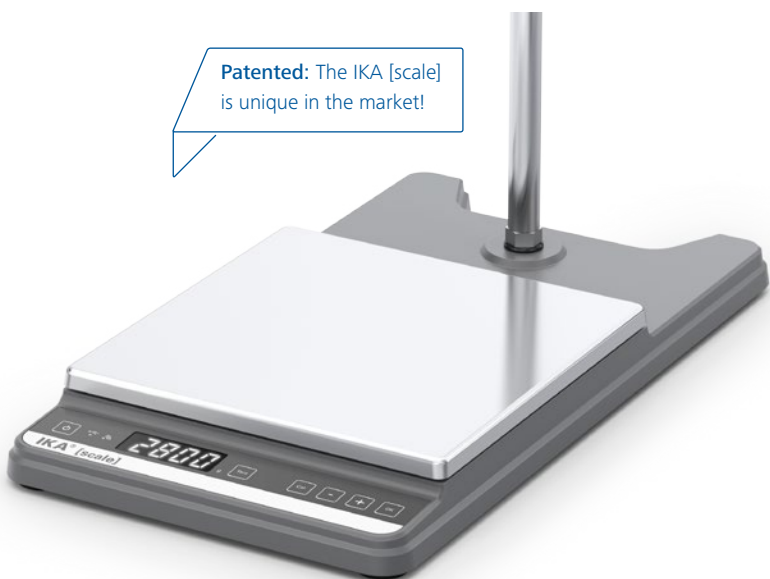

### Sample weight of up to 6 kg

The intelligent stand scale has a weighing range of up to 6 kg with a resolution of 1 g. It is therefore particularly suitable for the direct and controlled addition of sample material while stirring.

### Clear display

Working with the clear and modern LCD display is very user-friendly. All information, whether measurement values, calibration steps or status messages, is shown clearly.

### Stable set-up plate

The solid die cast zinc housing provides stability and robustness.

### Bluetooth connection

The weight of the sample measured can directly be transferred to a suitable device (e.g. IKA Mini-/Microstar series) and shown on its display.

#### TECHNICAL DATA

Overload protection	Yes
Dimensions (W × H × D)	215 × 45 × 390 mm
Weight	3.0 kg
Permissible ambient temperature	5 – 40 °C
Permissible relative humidity	80 %
Protection class according to DIN EN 60529	IP 21
Interface	USB
Voltage	100 – 240 V
Frequency	50/60 Hz

IKA [scale]

Ident. No. 0025004318

IKA-Werke GmbH & Co. KG

Janke & Kunkel-Straße 10, 79219 Staufen, Germany

Phone: +49 7633 831-0, Fax: +49 7633 831-98

eMail: sales@ika.de



www.ika.com



IKAworlwide // #lookattheblue